

Vacuum Pump V-700



The V-700 Vacuum Pump is a chemically resistant PTFE diaphragm pump designed for use with all gases and vapors typically associated in the lab. With a flow rate of 1.8 m³/h and an ultimate vacuum of below 10 mbar, it covers a wide range of classical laboratory applications and is optimally suited for work with a rotary evaporator. Other typical applications include the evacuation of drying ovens or vacuum centrifuges and the generation of vacuum for filtration purposes or for conveying liquids.


The vacuum pump offers the following advantages:

- Unique window head allows for visual inspection of diaphragm making it easy to check for solvent build-up or detection of contamination
- Innovative, single-stroke speed control for hysteresis-free and precise vacuum (in combination with the Vacuum Controller V-850/855)
- Sophisticated soundproofing with enclosed housing and counterbalancing minimize noise and vibration
- Integrated gas ballast for drying the diaphragms during operation, ensuring consistently high flow rates
- Virtually maintenance-free, easy to access diaphragms
- Compact design and solid workmanship, with integrated carrying handle
- Strictly axial movement of the diaphragm increases service life




Pump head design with window (blue) and diaphragm (green)

Minimum emission thanks to sophisticated secondary condenser, completely enclosed housing and Vacuum Controller

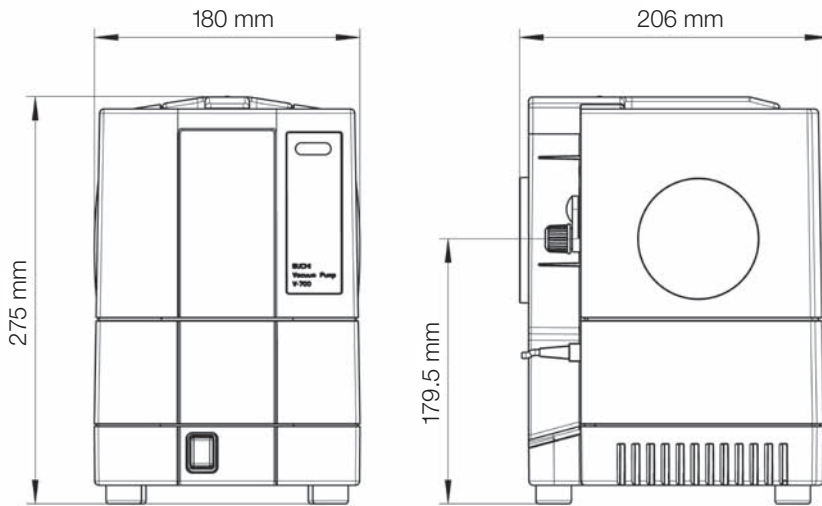


Volume flow rate (DIN 28432)	1.8 m ³ /h
Number of stages (heads)	2 (2)
Ultimate vacuum (absolute)	<10 mbar
Ultimate vacuum (with gas ballast)	24 mbar
Vacuum connection	GL 14
Power consumption	210 W
Supply voltage	100–240 V/50–60 Hz
Motor	brushless DC motor
Revolution speed	max. 1600 rpm min.-1
Sound level (DIN 45635)	40–52 dB (A) dependent on operation mode
Safety class	IP 34
Weight	5.3 kg
Materials in contact with medium	PEEK, PTFE, glass, FEP
Approvals	CE

Volume flow rate (m³/h) of a V-700 dependent on inlet pressure (mbar)

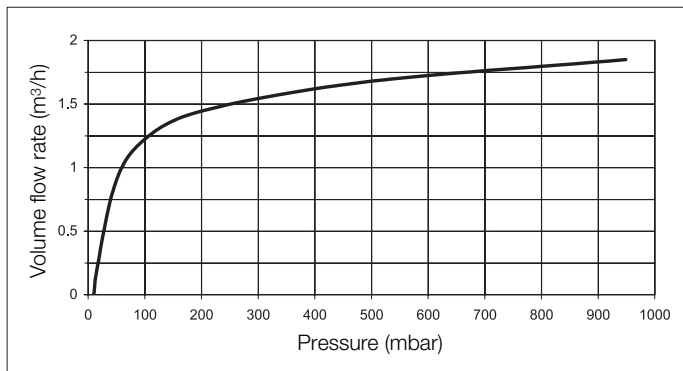

Evacuation time (s) of a Rotary Evaporator (Volume = 3.4 L) with a V-700

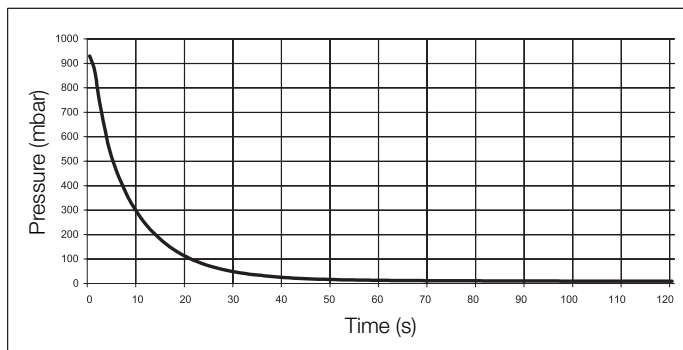

Ordering data:

V-700 Vacuum Pump	
85–264 V/50–60 Hz	
Order No.	71000
Scope of delivery:	
V-700 Vacuum Pump, sound absorber, main cable and operating instructions	

Accessories:

Accessories:	Order No.
Secondary condenser cpl.	47180
Secondary cold trap cpl.	47190
Vacuum Controller V-850 cpl. (incl. 47280)	47299
Vacuum Controller V-855 cpl. (incl. 47280)	47298
Manual vacuum regulation, cpl. (incl. 47280)	47291
Communication and fixation set	47280
Woulff' bottle cpl. incl. holder	47170
(Recommended in combination with V-850/855 or manual regulation)	
Cable RJ45 2000 mm	44989
Cable Mini-DIN 1500 mm	38010

BUCHI Labortechnik AG
Postfach
9230 Flawil 1
Schweiz
T +41 71 394 63 63
F +41 71 394 65 65
buchi@buchi.com
www.buchi.com

BUCHI Labortechnik GmbH
Postfach 10 03 51
45003 Essen
Deutschland
Freecall 0800 414 0414
T +49 201 747 490
F +49 201 237 082
deutschland@buchi.com
www.buechigmbh.de

BUCHI Labortechnik GmbH
Branch Office Netherlands
Postbus 142
3340 AC Hendrik-Ilo-Ambacht
The Netherlands
T +31 78 684 94 29
F +31 78 684 94 30
netherlands@buchi.com
www.buchi.nl

BUCHI Italia s.r.l.
Centro Direzionale, Milano Fiori
Pal. A-4, Strada 4
20090 Assago (MI)
Italia
T +39 02 824 50 11
F +39 02 57 51 28 55
italia@buchi.com
www.buchi.it

BUCHI (Thailand) Ltd.,
77/175, Sin Sathon Tower,
39th FL, Unit F
Krungthongburi Rd.
Klongtonnai, Klongsan
Bangkok 10600
Thailand
T +66 2 862 08 51
F +66 2 862 08 54
bacc@buchi.com
www.buchi.com

BUCHI SMP
Services Private Ltd.
201, Magnum Opus
Shantinagar Industrial Area
Vakola, Santacruz (East)
Mumbai 400 055,
India
T +91 22 66 98 94 50 / 51
F +91 22 66 98 94 52
smp@buchi.com
www.buchi.com

BUCHI Corporation
19 Lukens Drive, Suite 400
New Castle
Delaware 19720
USA
T +1 302 652 3000
F +1 302 652 8777
Toll Free: +1 877 692 8244
us-sales@buchi.com
www.mybuchi.com

BUCHI Hong Kong Ltd.
1810 Fortress Tower
250 King's Road
North Point, Hong Kong
China
T +852 2389 2772
F +852 2389 2774
china@buchi.com
www.buchi.com

BUCHI Shanghai Trading LLC
21/F Shanghai Industrial
Investment Building
18 Caoyi Bei Road
200030 Shanghai
China
T +86 21 6468 1888
F +86 21 6428 3890
china@buchi.com
www.buchi.com

BUCHI UK Ltd
5 Whitegate Business Centre
Jardine Way
Chadderton
Oldham OL9 9QL
United Kingdom
T +44 161 633 1000
F +44 161 633 1007
uk@buchi.com
www.buchi.co.uk

BUCHI Sarl
5, rue du Pont des Halles
Z.A. du Delta
94656 Rungis Cedex
France
T +33 1 56 70 62 50
F +33 1 46 86 00 31
france@buchi.com
www.buchi.fr

Nihon BUCHI K.K.
3F IMON Bldg.,
2-7-17 Ikenohata, Taito-ku,
Tokyo 110-0008
Japan
T +81 3 3821 4777
F +81 3 3821 4555
nihon@buchi.com
www.nihon-buchi.co.jp

We are represented by more than 100 distribution partners worldwide. Find your local representative at

www.buchi.com

Quality in your hands